DCOØK

BÁSCULA MIDE GRASA
BALANÇA MEDE GORDURA
BODY FAT SCALES
PÈSE-PERSONNE IMPÉDANCEMÈTRE

MANUAL DE INSTRUCCIONES
MANUAL DE INSTRUÇÕES • INSTRUCTION MANUAL
NOTICE D'UTILISATION

BÁSCULA MIDE GRASA

DESCRIPCIÓN

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
BODY FAT SCALES

DESCRIPTION

1. Weighing platform
2. Weight selector button
3. Battery compartment
4. LCD display
5. Measuring sensors
6. Person selector: age, sex, and height

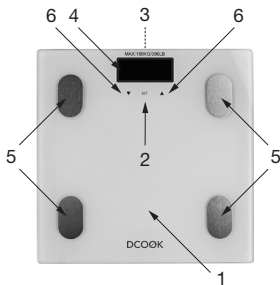SAFETY INSTRUCTIONS

• **Read these instructions carefully before switching on the appliance and keep them for future reference. Failure to follow and observe these instructions could lead to an accident.**

• **IMPORTANT:** This appliance can be used by children over 8 years and people with reduced physical, sensorial or mental capacities or with a lack of experience or knowledge if they are supervised or have received training on the safe use of the device and the risks involved.

- Children should not play with the appliance and must be supervised.
- Keep the appliance out of the reach of children under the age of 8 years.
- Check the appliance for damage from time to time. Never use the appliance if the glass is cracked or broken. In such a case, repairs must be carried out by a qualified electrician (*). Avoid any blows and falls.
- Take care not to exceed the maximum weighing capacity, as, otherwise, it could be damaged.
- Never use the appliance outdoors, or subject it to weather conditions like rain, snow, etc. which could affect its proper working.
- Never immerse the appliance in water or any other liquid, as it could stop working.
- Do not step on the appliance with wet feet and keep the surface dry.
- This is an appliance for domestic use. It is not suitable for industrial or commercial use.

Competent qualified electrician: Any person who is qualified, authorised and competent to do this type of repairs with an aim to avoid danger of any kind.

WARNING CONCERNING THE BATTERIS

- Change the batteries when they are flat.
- Only use the batteries of the same model or one equivalent to it.
- Remove the battery from the compartment if you do not intend to use the appliance over a long period of time.
- The batteries must be placed in the special containers provided for them.
- Remove the batteries from the appliance if you wish to dispose of the scale.
- If you detect any kind of leak in the batteries, do not handle them with your bare hands. Cover them.
- Always place the batteries in their compartment, with the correct polarity.
- Never throw the batteries into the fire or expose the appliance to the sun.
- Never mix new batteries with used or old ones.
- The battery compartment must always be far from fire or any other source of heat.

HOW TO USE

1. Step onto the scale for a simple weight measurement.
2. When the scale is turned on, it can measure weight, body fat, water content, muscle mass, bone mass, calories and BMI.
3. Touch the setting button so that the display window will enter the current user data. At this time, you can measure weight, body fat, water content, muscle mass, bone mass, calories, and BMI on the scale. You can also reset the user data by touching the setting button.
4. You need to set height, sex and age manually. These are stored after being set. The next time you use the scale directly.
 5. This scale can store up to 10 user profiles.
 6. The scale turns off automatically when not in use.

Data Lock function: the data will stay on the screen for a few seconds after stepping off the scale.

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1. User number

Turn on the scale. When the scale displays 0.00 kg, press “ SET ”. The user number symbol will flash, use “ ▲ ” or “ ▼ ” to select the user number and press “ SET ” to confirm.

2. Height

After confirming the user number, the height value flashes. Use “ ▲ ” or “ ▼ ” to set the height value and press “ SET ” to confirm.

3. Age

After confirming the height, the age value flashes. Use “ ▲ ” or “ ▼ ” to set the age and press “ SET ” to confirm.

4. Sex

After confirming the age, the sex symbol flashes. Use “ ▲ ” or “ ▼ ” to select the sex and press “ ” to confirm. After the user parameters are set, the scale shows 0.00 kg. Now you can measure weight, body fat, water content, muscle mass, bone mass, calories, and BMI. Stand barefoot on the scale with your feet in contact with the steel part. Remain still. The scale will display your weight, body fat, water content, muscle mass, bone mass, calories, and BMI one by one.

Body fat level reference chart

		Body fat percentage (%)			
	Age	Slightly thin	Average	Slightly fat	Fat
Man	≤ 25	< 10%	10,1 - 21%	21,1 - 26%	> 26,1%
	> 25	< 13%	13,1 - 23%	23,1 - 28%	> 28,1%
Woman	≤ 25	< 20%	20,1 - 34%	34,1 - 39%	> 39,1%
	> 25	< 21%	21,1 - 35%	35,1 - 40%	> 40,1%

Water content level reference chart

		Water content range		
	Age	Low	Average	Hight
Man	≤ 25	< 49%	49 - 70%	> 70%
	> 25	< 47%	47 - 68%	> 68%
Woman	≤ 25	< 50%	50 - 75%	> 75%
	> 25	< 45%	45 - 70%	> 70%

Muscle mass level reference chart

		Muscle mass range		
	Age	Low	Average	High
Man	≤ 25	< 37%	37 - 64%	> 64%
	> 25	< 35%	35 - 59%	>59%
Woman	≤ 25	< 33%	33 - 55%	> 55%
	> 25	< 30%	30 - 50%	> 50%

Bone mass level reference chart

		Bone mass range		
	Age	Low	Average	High
Man	≤ 25	< 3%	3 - 9%	> 9%
	> 25	< 3%	3 - 8,5%	> 8,5%
Woman	≤ 25	< 2,5%	2,5 - 8%	> 8%
	> 25	< 2,5%	2,5 - 7%	>7%

Body mass index (BMI)

Slightly thin	Average	Slightly fat	Fat
< 18,5	18,5 - 24,9	25 - 29	> 30

CARE AND CLEANING:

- This appliance must be cleaned with a moist cloth and a few drops of detergent.
- Do not use abrasive or sharp products to clean this appliance.
- Never immerse this appliance in water or any other liquid.

- Do not store the scale in a vertical position. Always keep it in a horizontal position.

TECHNICAL SPECIFICATIONS

- Model: 5430172
- 3V DC 3V(2AAAx1.5V)
- 180kg/100gr (396Lb/0.2Lb)
- Muscle mass measurement range: 25-75%
- Height range: 100-240 cm
- Bone mass measurement range: 0.1-10%
- Age range: 10-99
- Water content measurement range: 20-80%
- High-precision sensors
- Body fat measurement range: 5-60%
- Automatic step-on and shut-off
- Calories measurement range: 500-5000
- Body fat percentage error display: "Err"
- Low battery display: "Lo"
- Overweight display: "oL"

This appliance conforms to Directive 2014/30/EU on Electromagnetic Compatibility.

ENVIRONMENTAL PROTECTION - / EU DIRECTIVE 2012/19

In order to preserve our environment and protect our health, any waste electrical and electronic equipment must be disposed of in accordance with specific rules with the implication of suppliers and users.

The symbol on the packaging indicates that this product must not be disposed of together with ordinary household waste. The user must take the product to a collection point for the recycling of electrical and electronic appliances.

Disposal must be in line with prevailing local regulations. Consult your regional or municipal authorities about the options for waste disposal.

GUARANTEE TERMS

The legal guarantee conforming to prevailing legislation is applicable to this product. For the duration of the legal guarantee, all material and manufacturing defects are covered by the guarantee.

To avail of the guarantee, you will need to present the purchase receipt at the store where the item was bought. You can then choose either a refund or an exchange for a new item.

For after-sales service contact us with any query you may have at: <http://adibera.com/es/contact/questions>

No claims will be accepted for broken accessories or defects due to inappropriate use or care, excessive wear due to intensive, non-domestic use, incorrect use due to failure to comply with the instruction manual, or normal wear. If this appliance is used incorrectly, for professional or semi-professional purposes or in a way that is not in line with the operating instructions, the guarantee will be null and void.

No guarantee claims will be accepted if the appliance has been modified or if any attempts have been made to repair it outside the official technical assistance services. In such an event, all liability in the event of accidents will be rejected.

The parts subjected to wear, as well as cleaning and maintenance or the replacement of worn parts are not covered by the guarantee and must be paid for by the buyer.

This instruction manual and any future updates can be digitally downloaded at www.adibera.com

AFTER THE LEGAL PERIOD OF GUARANTEE

Repairs after the guarantee period will be on the owner's account.

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