# DCOOK

**BÁSCULA MIDE GRASA**  
BALANÇA MEDE GORDURA  
BODY FAT SCALES  
PÈSE-PERSONNE IMPÉDANCEMÈTRE

**MANUAL DE INSTRUCCIONES**  
MANUAL DE INSTRUÇÕES • INSTRUCTION MANUAL  
NOTICE D'UTILISATION

# BÁSCULA MIDE GRASA

## DESCRIPCIÓN

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
# BODY FAT SCALES

## DESCRIPTION

1. Weighing platform
2. Weight selector button
3. Battery compartment
4. LCD display
5. Measuring sensors
6. Person selector:  
age, sex, and height
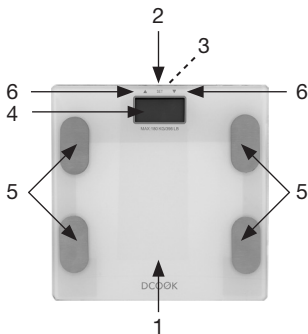
## SAFETY INSTRUCTIONS

• **Read these instructions carefully before switching on the appliance and keep them for future reference. Failure to follow and observe these instructions could lead to an accident.**

• **IMPORTANT:** This appliance can be used by children over 8 years and people with reduced physical, sensorial or mental capacities or with a lack of experience or knowledge if they are supervised or have received training on the safe use of the device and the risks involved.

- Children should not play with the appliance and must be supervised.
- Keep the appliance out of the reach of children under the age of 8 years.
- Check the appliance for damage from time to time. Never use the appliance if the glass is cracked or broken. In such a case, repairs must be carried out by a qualified electrician (\*). Avoid any blows and falls.
- Take care not to exceed the maximum weighing capacity, as, otherwise, it could be damaged.
- Never use the appliance outdoors, or subject it to weather conditions like rain, snow, etc. which could affect its proper working.
- Never immerse the appliance in water or any other liquid, as it could stop working.
- Do not step on the appliance with wet feet and keep the surface dry.
- This is an appliance for domestic use. It is not suitable for industrial or commercial use.

Competent qualified electrician: Any person who is qualified, authorised and competent to do this type of repairs with an aim to avoid danger of any kind.

## WARNING CONCERNING THE BATTERIS

- Change the batteries when they are flat.
- Only use the batteries of the same model or one equivalent to it.
- Remove the battery from the compartment if you do not intend to use the appliance over a long period of time.
- The batteries must be placed in the special containers provided for them.
- Remove the batteries from the appliance if you wish to dispose of the scale.
- If you detect any kind of leak in the batteries, do not handle them with your bare hands. Cover them.
- Always place the batteries in their compartment, with the correct polarity.
- Never throw the batteries into the fire or expose the appliance to the sun.
- Never mix new batteries with used or old ones.
- The battery compartment must always be far from fire or any other source of heat.

## HOW TO USE

- Before using the appliance, remove the protective plastic from the batteries in order to place them in the scale to get it working.
- Exert slight pressure on the weight selector to choose the required weight unit, Kg or Lb.
- The selected symbol, Kg or Lb, will display on the LCD.
- To get the correct results, you are advised to weigh yourself with no shoes and make sure your feet are in contact with the 4 metallic sensors on the weighing platform.
- This appliance has 10 memories, for 10 different users.
- Each user is differentiated by the letters P0 to P9.
- To memorise each user, take the following steps:

- 1- Press SET, and P0...P9 will display on the screen. They will flash until you select the required one.
- 2- Use the arrows to select the required memory and press SET again to store it and move on to the next step where you must select the sex (male or female) according to the following silhouettes:

- 3- Use the arrows to select and store the sex and move on to the next step and press SET
- 4- You must now select the age (AGE), by using the arrows. Once selected, press SET and move on to set the height. Press SET again to store the height.

- Once all the data has been stored in the memory, its use is as simple as with any ordinary scale.
- Press SET and select the memory with the arrows. Once the scale displays 0.0, place the appliance on a flat, smooth and horizontal surface. Step on the appliance with dry, bare feet, making sure to cover the 4 sensors and stay still without touching anything else.

The appliance will immediately display the following information intermittently:

1- The weight in kg/lb

2- % of body fat 

3- Body moisture 

4- Muscle mass 

5- Bone proportion 

- The results will show on the display twice and then it will disconnect.
- **Note:** the information shown in the table below is for indicative purposes only. If you have any queries about the proportion of body fat, please consult your doctor.

|       | Age     | Low    | Recommended | High    | Very high |
|-------|---------|--------|-------------|---------|-----------|
| Woman | 20 - 39 | 5 - 20 | 21 - 33     | 34 - 38 | >38       |
|       | 40 - 59 | 5 - 22 | 23 - 34     | 35 - 40 | >40       |
|       | 60 - 79 | 5 - 23 | 24 - 36     | 37 - 41 | >41       |
| Man   | 20 - 39 | 5 - 7  | 8 - 20      | 21 - 25 | >25       |
|       | 40 - 59 | 5 - 10 | 11 - 21     | 22 - 27 | >27       |
|       | 60 - 79 | 5 - 12 | 13 - 25     | 26 - 30 | >30       |

- **The auto power off function:** following a few seconds of inactivity, the appliance will switch itself off automatically, to save energy.
- **LO message:** if this message displays on the screen, it means that the batteries need to be changed soon.
- **ERR message:** this message means that the maximum weighing capacity has been exceeded.

## CARE AND CLEANING:

- This appliance must be cleaned with a moist cloth and a few drops of detergent.
- Do not use abrasive or sharp products to clean this appliance.
- Never immerse this appliance in water or any other liquid.
- Do not store the scale in a vertical position. Always keep it in a horizontal position.

## TECHNICAL SPECIFICATIONS

Model: 8084001

3V DC ( 1 x CR2032)

180kg/100gr ( 396Lb/0.2Lb)

**This appliance conforms to Directive 2014/30/EU on Electromagnetic Compatibility.**

## ENVIRONMENTAL PROTECTION - / EU DIRECTIVE 2012/19

In order to preserve our environment and protect our health, any waste electrical and electronic equipment must be disposed of in accordance with specific rules with the implication of suppliers and users.

The symbol on the packaging indicates that this product must not be disposed of together with ordinary household waste. The user must take the product to a collection point for the recycling of electrical and electronic appliances.

Disposal must be in line with prevailing local regulations. Consult your regional or municipal authorities about the options for waste disposal.

### GUARANTEE TERMS

The legal guarantee conforming to prevailing legislation is applicable to this product. For the duration of the legal guarantee, all material and manufacturing defects are covered by the guarantee.

To avail of the guarantee, you will need to present the purchase receipt at the store where the item was bought. You can then choose either a refund or an exchange for a new item.

For after-sales service contact us with any query you may have at: <http://adibera.com/es/contact/questions>

No claims will be accepted for broken accessories or defects due to inappropriate use or care, excessive wear due to intensive, non-domestic use, incorrect use due to failure to comply with the instruction manual, or normal wear. If this appliance is used incorrectly, for professional or semi-professional purposes or in a way that is not in line with the operating instructions, the guarantee will be null and void.

No guarantee claims will be accepted if the appliance has been modified or if any attempts have been made to repair it outside the official technical assistance services. In such an event, all liability in the event of accidents will be rejected.

The parts subjected to wear, as well as cleaning and maintenance or the replacement of worn parts are not covered by the guarantee and

must be paid for by the buyer.

This instruction manual and any future updates can be digitally downloaded at [www.adibera.com](http://www.adibera.com)

## **AFTER THE LEGAL PERIOD OF GUARANTEE**

Repairs after the guarantee period will be on the owner's account.

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